



ViscoDispo

The ViscoDispo disposable sampler can be used to take all layer and target point samples of high-viscosity media up to 100.000 mPas. The ergonomically shaped handle and the very stable drawbar make it easy to draw even viscous substances.

- Ergonomic handle
- No seals
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-3320	500 mm	100 mL	standard	20
853933-3321	500	100	sterilized	20
853933-3330	1000	200	standard	20
853933-3331	1000	200	sterilized	20



LiquiDispo

The LiquiDispo disposable sampler can be used to take all layer and target point samples sample of liquids from highly liquid to slightly viscous. An ergonomic handle makes it easier to draw the sample: up to 200 ml can be drawn, as with a syringe.

- Ergonomic handle
- No seals
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-1120	500 mm	100 mL	standard	20
853933-1121	500	100	sterilized	20
853933-1130	1000	200	standard	20
853933-1131	1000	200	sterilized	20

www.dynalon.com



DispoTube

The ergonomic handle of the DispoTube makes it easier to use and facilitates a good grip even under difficult conditions, for example with gloves. Disposable sampler for taking samples of liquids.

- Ergonomic handle
- All-layer and target point samples possible
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-6620	500 mm	100 mL	standard	20
853933-6621	500	100	sterilized	20
853933-6630	1000	200	standard	20
853933-6631	1000	200	sterilized	20



DispoPipette

Simple disposable sampler for taking samples of low viscosity media. If an opened DispoPipette is inserted into the liquid an all-layer sample is obtained. If it is closed while it is inserted up to the sampling point and then opened, a targetpoint sample is obtained. Very easy to operate. With two opening sizes for different viscosities.


- All-layer and target point samples possible
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-5500	500 mm	100 mL	standard	20
853933-5501	500	100	sterilized	20
853933-5531	1000	200	standard	20
853933-5532	1000	200	sterilized	20

LaboPlast® | SteriPlast® DispoSampler

NEW

- ▶ For disposable use
- ▶ Individually packaged
- ▶ Production, assembly and packaging in class 7 clean room (10,000)
- ▶ With EU foodstuffs and FDA approval
- ▶ Optionally sterilised by gamma rays

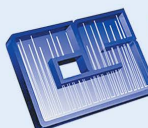

**CHROMATOGRAPHIC
SPECIALTIES INC.**



1-800-267-8103 • www.chromspec.com • tech@chromspec.com



DispoSampler – disposable sampler | for the highest standards in sampling

Disposable samplers – DispoSampler – are used to take samples of liquids, viscous media, powders and granulates. Target point, all-layer and cross-section samples can be taken, as well as small and large sample quantities. Unlike conventional reusable samplers, disposable samplers do not need to be cleaned after use.

- Cleaning reusable samplers is a demanding and often complicated task, with greasy, strongly adhesive or toxic substances. Sometimes the task is even impossible.
- Work processes do not have to be interrupted by cleaning
- Cross contamination is completely eliminated
- The cleaning method does not need to be validated in audits
- Time and money are saved

www.dynalon.com



MultiDispo

Disposable zone samplers for taking samples from multiple layers of the sample material. Suitable for free-flowing powders and small-grain granulates. The representative all-layer sample can be emptied directly through the open handle into a sample container.

- Easy emptying through the opening in the handle
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-4440	500 mm	100 mL	standard	20
853933-4441	500	100	sterilized	20
853933-4450	1000	200	standard	20
853933-4451	1000	200	sterilized	20



ZoneDispo

The ZoneDispo disposable all-layer sampler is used in sampling to collect a representative all-layer sample from three different layers. Suitable for free-flowing powders and small-grain granulates. The collected sample can be secured and transported in the closed sampler.

- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-4420	500 mm	100 mL	standard	20
853933-4421	500	100	sterilized	20
853933-4430	1000	200	standard	20
853933-4431	1000	200	sterilized	20



TargetDispo

TargetDispo is used to take a target point sample at a position in the sample material that has previously been precisely determined. The sampler pierces into the powder or granulate while closed. When it reaches the sample removal point it is opened and then closed again. In this way a precise target-point sample is obtained. The sample can simply be emptied through the open handle.

- Easy emptying through the opening in the handle
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-4460	500 mm	100 mL	standard	20
853933-4461	500	100	sterilized	20
853933-4470	1000	200	standard	20
853933-4471	1000	200	sterilized	20

PowderDispo

PowderDispo is the ideal disposable sampler for cross-section samples of powdery or slightly cohesive substances. A complete all-layer sample can be taken from bags or smaller containers with the lance.

- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-5510	300 mm	60 mL	standard	20
853933-5511	300	60	sterilized	20
853933-5522	500	100	standard	20
853933-5521	500	100	sterilized	20



DispoPicker

Disposable sampler for powders and small-grain bulk materials. Because the tip is closed, the sample material cannot trickle back out of the sample lance. DispoPicker is ideally suited for cross-section samples.

- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-5540	300 mm	60 mL	standard	20
853933-5541	300	60	sterilized	20
853933-5550	500	100	standard	20
853933-5551	500	100	sterilized	20



DispoLance

The DispoLance disposable sampler can be used for taking large and small sample quantities quickly and easily, directly from bags. The sampler is suitable for free-flowing powders and small-grain granulates. The DispoLance sampler pierces into the container with the sturdy tip and the sample flows through the open inner tube directly into the container or bag that is held underneath.

- Variable sample quantity
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Model	Pack Quantity
853933-5560	300 mm	standard	20
853933-5561	300	sterilized	20
853933-5570	500	standard	20
853933-5571	500	sterilized	20



MicroDispo

Especially small sample quantities (10 ml) can be taken with MicroDispo. In pharmaceutical applications especially, often only a small quantity can or should be extracted for analysis. Using the disposable sampler, precise target point samples can be taken of free-flowing powders and small-grain bulk materials.

- Ergonomic handle
- HDPE
- Clean room manufactured
- Optionally sterilized by gamma rays

Item	Length	Volume	Model	Pack Quantity
853933-2240	500 mm	10 mL	standard	20
853933-2241	500	10	sterilized	20
853933-2250	1000	10	standard	20
853933-2251	1000	10	sterilized	20

