
**RESTEK INNOVATION** 

### **NEW Restek Electronic Leak Detector**



Why allow a tiny leak become a GIANT

problem?

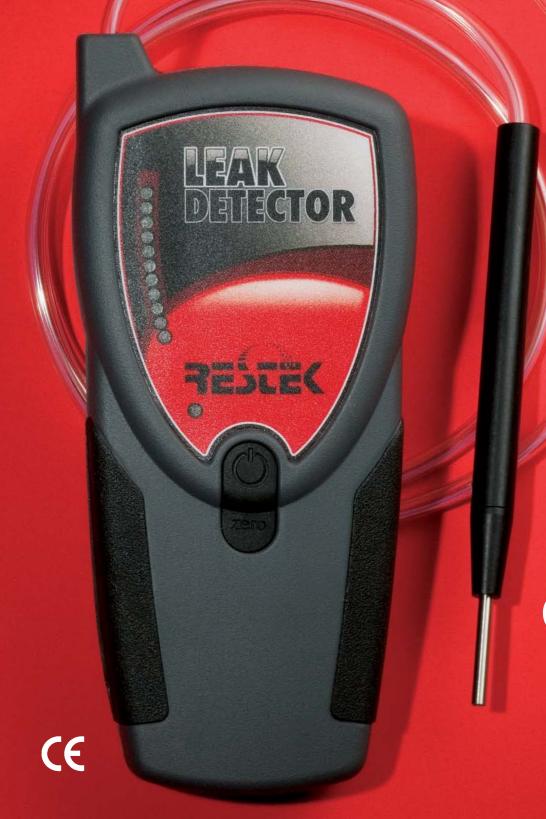

**Chromatography Products** 



# Contact your local Restek representative to **order yours TODAY!**





Carrying/storage case included with purchase of unit.

Avoid using liquid leak detectors on a capillary system! Liquids can be drawn into the system.

Caution: The Restek Electronic Leak Detector is NOT designed for determining leaks in a combustible environment. A combustible gas detector should be used for determining combustible gas leaks under any condition. The Restek Electronic Leak Detector may be used for determining trace amounts of hydrogen in a GC environment only.

## NEW Restek Electronic Leak Detector

Ever have a small leak turn into a costly repair?

Protect your data and analytical column by using a Restek Leak Detector!

- Optimized sample flow path.
- A sleek, new ergonomic, hand-held design.
- Rugged side grips for added durability.
- Handy probe storage for cleanliness.
- Longer lasting battery, up to 6 hours of continuous use.
- Automatic shut-off capabilities.
- A convenient carrying and storage case.
- A universal charger set (US, European, UK and Australian plugs included).

Backed by a 1 year warranty, the new Restek Leak Detector will again set an industry standard for performance and affordability in a hand-held leak detector.

#### **Leak Detector Facts**

Ecun Detector ructs	
Detectable Gases:	helium, nitrogen, argon, carbon dioxide, hydrogen
Battery:	Rechargeable Ni-MH internal battery pack (6 hours normal operation)
Operating Temp. Range:	32°-120°F (0°-48°C)
Humidity Range:	0-97%
Warranty:	one year
Certifications:	CE, Japan
Compliance:	WEEE, RoHS

#### **Detection Limits**

	Minimum Detectabe Leak	Indicating LED Light Color	
Gas	Rate (atm cc/sec.)		
Helium	1.0 X 10 <sup>-5</sup>	Red	
Hydrogen**	1.0 X 10 <sup>-5</sup>	Red	
Nitrogen	1.4 X 10 <sup>-3</sup>	Yellow	
Argon	1.0 X 10 <sup>-4</sup>	Yellow	
Carbon Dioxide	1.0 X 10 <sup>-4</sup>	Yellow	

### Restek's Competitive Advantage

	Restek Leak Detector	Agilent 63388A Leak Detector
Probe Storage	<b>✓</b>	
Ergonomic Design	<b>✓</b>	
Hand Grips	<b>✓</b>	
Optimized Sample Flow Path	<b>✓</b>	
Storage Case	<b>✓</b>	
LED Leak Strength Indicator	V	
6-Hour Battery	<b>✓</b>	



Soft-side storage case is ideal for storing your leak detector in smaller spaces such as your tool box.



Verify hard-to-reach leaks using the small probe adaptor.

Description	qty.	cat.#	
Leak Detector with Universal Charger Set (US, UK, European,			
Australian)	ea.	22839	
Soft-Side Carry/Storage Case	ea.	22657	
Small Probe Adaptor	ea.	22658	



Lit. Cat.# GNPC1078-INT

© 2008 Restek Corporation.

