


from Metronics®

## Dual Layer Color-Coded 1/16" PEEK® Tubing

- Inner layer of natural PEEK, outer layer of color-coded PEEK
- No possibility for stream to contact coloring agents
- Pressure performance equivalent to single layer tubing
- IDs from .004" to .030", ±.001"

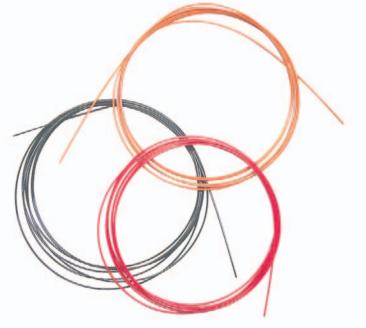

## **General Description**

Our Dual Layer PEEK tubing eliminates any concern that a critical sample stream could be contaminated by pigments used to color code the tubing. It looks like any other color-coded tubing at first glance, but a closer look reveals that the pigmented layer surrounds a separate but integrally-bonded inner layer of natural PEEK.

## **Specifications**

The 1/16" OD tubing is available in OEM quantities with any ID you specify, from .004" to .030". Standard sizes are listed below. The ID tolerance of ±.001" is equivalent to or better than that of tubing typically sold as premium grade. The tubing is available in a wide range of colors, striped or solid.\*

## 1/16" OD Standard Sizes

Color
Red
Yellow
Blue
Orange
Green



<sup>\*</sup> All colorants used in the manufacture of our tubing are RoHS-compliant (Reduction of Hazardous Substances)